VOLITO

(ES) Andador "VOLITO" - RF-1301

5 - 14

(IT) **(MT)** Deambulatore VOLITO - RF-1301

15 - 24

(PT) Andarilho "VOLITO" com Rodas - RF-1301

25 - 34

(GB) **(MT)** "VOLITO" Rollator - RF-1301

35 - 44

(DE) **(AT)** **(CH)** Rollator VOLITO - RF-1301

45 - 54

Este producto es conforme a la directiva 93/42/EWG de productos médicos.
Questo prodotto è conforme alla norma 93/42EWG sui prodotti medicinali.
Este produto está em conformidade com a Directiva 93/42/CEE sobre produtos de medicina.
This product conforms to the Medicinal Product guideline 93/42/EWG.
Dieses Produkt ist konform mit der Richtlinie 93/42/EWG über Medizinprodukte.

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ISTRUZIONI D'USO

ES **Manual De Instrucciones**

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pt@nws-service.com

www.nws-service.com

Read these instructions thoroughly. Retain for future reference! When passing this article on, please include these instructions.

Included in Delivery / Description of Parts

- 1x Rollator
- 1x Bell
- 1x Tray
- 1x Back belt with screw 3 x 10 mm
- 1x Rain hood
- 1x Upper stick holder
- 1x Lower stick holder
- 1x Instructions

fig. 1

- Ⓐ Hand grip
- Ⓑ Back belt with screw 3 x 10 mm
- Ⓒ Rain hood
- Ⓓ Upper stick holder
- Ⓔ Basket
- Ⓕ Lower stick holder
- Ⓖ Kerb deflector
- Ⓗ Tipping aid
- ① Securing screws
- Ⓜ Tray
- Ⓚ Brake lever
- Ⓛ Bell

Technical Details

- Max. user weight: 150 kg
- Max. basket load: 5 kg
- Max. tray load: 5 kg
- Max. length: 675 mm
- Max. width: 635 mm
- Rollator height: from 790 to 990 mm
- Folded dimensions: H 825 x W 280 x L 655 mm
- Rollator weight: approx. 7.8 kg (not including accessories)

- Hand grip length: 105 mm
- Hand grip width: 53 mm
- Distance to brake lever handle: 72 mm
- Tyre width: 33 mm
- Wheel diameter: 200 mm
- Seating height: 650 mm
- Seating width: 400 mm
- Turning circle: 860 mm

Symbols used

Observe the instructions

Date of manufacture

Batch designation

Manufacturer

Max. user weight

Max. tray load

Max. basket load

Max. width

Push handles height adjustment

For use outdoors and in enclosed rooms

Intended Use

The Rollator is to be used exclusively as a walking aid to promote the freedom of movement of physically affected people. It is suitable for outdoor use as well as in enclosed rooms. This article is conceived for private use and is unsuitable for commercial application.

Keep the Rollator close in to the body, whilst maintaining an upright position, making sure to always look forwards.

Never run with the Rollator and do not push the article too far ahead in order to avoid falling over or other dangers due to adverse weight.

Safety notes

Danger of injury! **WARNING!**

Danger of strangulation for children!

Never allow children to play unsupervised with this article or its packaging.

- Keep children away from this article. They could attempt to clamber onto it or pull themselves up on it, thereby causing injury.
- Use the article only for its stated use.
- The article should only be used by persons aware of these safety instructions.
- This article should not be used by persons exceeding the maximum user weight of 150 kg.
- Check the article before every use for damage or signs of wear and tear.
- Check that all of the screws are tight before every use!
- This article should only be used in whilst in good order and condition!
Do not use the article if you can see that it has been damaged or if you think it may be damaged!
- Use this article only when surface temperatures are suitable for its use.
Hot water, direct sunlight, sauna conditions, extreme cold, etc. can cause the surface of the article to assume temperatures which could endanger health.

- When parking the Rollator, make sure to use the brake so that it does not unintentionally roll away.
- If you are not using the Rollator, make sure that the brake is on so that no-one can fall over it, or that it does not unintentionally roll away.
- Only use the Rollator on a solid, flat surface. The tyres can catch in grooves and uneven recesses in the ground. Avoid drains, manhole covers and similar obstacles.
- Never stand on the Rollator and never attempt to carry passengers.
- Only use the Rollator when it has been completely and correctly assembled! Failure to observe this can cause injury to the user, and to others.
- The Rollator should only be used when in a completely unfolded condition.

Assembly / Settings

WARNING!


Only use the Rollator when it has been completely and correctly assembled! It should only be used in whilst in good order and condition! Do not use the article if you can see that it has been damaged or if you think it may be damaged!

Assembly

1. Remove the packaging materials and place the Rollator on the ground on its wheels. (fig. 2).

2. Press the seat down until it slots in with a click <CLICK> (fig.3, 4).

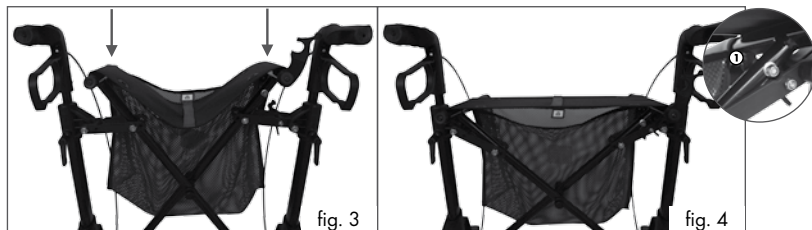
Bell

Note:

The bell can be mounted on either the right or left side.

Fasten the bell so that the ratchet ring passes through the intended opening. Remove the bell by pressing the release lever ① and pulling the ratchet ring out (fig.5, 6).

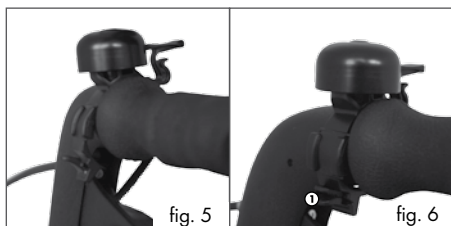Stick holder

Push the stick holder up ② and down ③ in the manner intended (fig.7).

Setting the hand grip height


⚠ Warning!

- **Make sure that the push handles are set at the same height in order to maintain safe transport and use.**
- **The push handles must click audibly into the 7 holes ③ <CLICK>.**
- **Turn the retaining screws only up to the red line.**

Note:

We recommend the hand grip height to be set at the same height as that of the wrists when hanging down at arm's length when standing.

1. In order to set these at the correct height, release the push handle retaining screws ①, (fig.8).
2. Pull the push handles out to the required height and then tighten the retaining screws. If the retaining screws are not exactly parallel with the frame tubing after tightening, press the locking button in ② and arrange the retaining screws in the correct position (Fig.9).

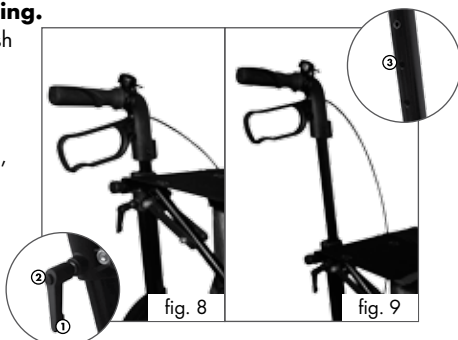Brakes

WARNING!

Always use both brakes simultaneously. Braking effect can be hindered by dirty, wet, icy, slippery or snowy surfaces.

Note:

Do not attempt to push the Rollator whilst the brakes are activated and do not pull the brakes over a lengthy period of time when pushing so as to avoid excessive wear and tear on the tyres.

• **Braking operation**


To brake when on the move, pull the brake lever upwards, towards the handle (fig.10).

• **Locking brake**

Apply the locking brake when the Rollator is not in use. This will ensure that it does not roll away unintentionally.

1. Press the brake lever downwards, over the point at which pressure can be felt.
The Rollator is now protected from rolling away.
2. To release the locking brake, pull the brake lever upwards, over the point at which pressure can be felt, until the brake has been released (fig.11).

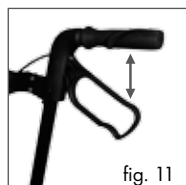Back belt

Note

For assembly and disassembly you will require a cross-head screwdriver. The 3 x 10mm screws are slightly turned into their respective holes in delivered default condition.

1. Push the back belt with the 3 x 10mm screws from both sides and from above into the slots ① in the hand grips (fig.12).

2. Tighten the 3 x 10mm screws on both sides in order to secure the back belt (figs. 13, 14).

3. To dismantle the back belt, loosen the 3 x 10mm screws from their holes and pull the back belt out of the slots in an upward direction. Re-tighten the screws into their holes in the back belt so that they do not get lost.

Basket

WARNING!

- The basket is suitable only for transporting small objects.
- The maximum load should not exceed 5kg.

1. Assembling/Dismantling the basket is achieved with the use of 2 push-buttons and hook and loop fasteners on both sides of the basket (fig.15).

2. Assemble the rain hood by attaching it to the push-buttons on both sides of the basket (Fig.16).

3. If necessary, pull the rain hood over the basket using the straps ② (fig.17, 18).

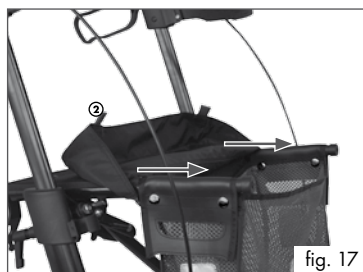
fig. 17

fig. 18

Seat

⚠ WARNING!

- **The Rollator seat is not intended to support weight over a lengthy period of time and is not intended for use as a chair.**
- **The seat is currently not intended for more than one person.**

If you intend to use the seat, please observe the following before you sit down:

- **The Rollator must be standing on a solid, flat surface; otherwise there is a risk of it rolling away.**
- **Always activate both locking brakes before you sit down.**
- **Sit down in the Rollator so that you face away from the direction of travel. You can lie back against against the back belt.**
- **Always sit upright in the seat. Your feet should always remain in contact with the ground.**

Tray

⚠ WARNING!

- **The tray is not heat-resistant and should not be used to carry heavy, sharp or fragile objects.**
- **The maximum load the tray can carry is 5kg. Spread the load evenly over the tray.**
- **Do not sit on the seat when the tray is mounted!**

1. Use the tray supplied and place it on the seat of the Rollator.

The pins on the underside of the tray must glide noticeably into the recesses ① of the seating area (fig.19).

fig. 19

Folding the Rollator


WARNING!


Danger of crushing!

When folding the Rollator, make sure your hands do not get caught in the moving parts. This can cause injury.


1. Remove the tray from the Rollator.
2. Unlock the fold lock ① by pulling the lever under the seat upwards. (fig.20).

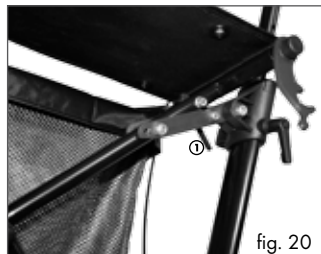3. Hold on to the Rollator's seat loop ② and pull it upwards (fig.21).

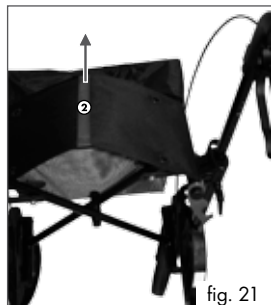4. Then secure the Rollator with the locking clamp ③ (fig.22, 23).

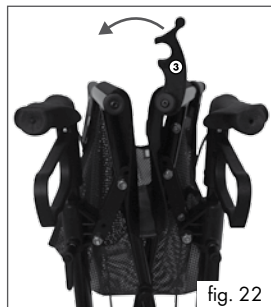Care and Maintenance

Care

CAUTION!

Do not use any abrasive or corrosive cleaning agents or sponges. These can damage the surface. Never use a high-pressure cleaner! This can also damage the surface of the Rollator.

1. Wipe over the Rollator at regular intervals with a damp cloth and perhaps a mild detergent when necessary.
2. Make sure that you dry it thoroughly, especially the handles before use or storage.


Maintenance

WARNING!

Do not use the Rollator when the tyres, brakes or frame appear to be damaged!

Brakes

1. Check the brakes before each use to make sure they still work. If they are not working, do not attempt to repair them yourself. Contact your dealer or service centre.
2. Check the brake cables for tears or kinks.
3. Calibrate: Set the distance between the brake pads and the wheel to approximately 2mm. Release the counter nut ② and adjust the brake cable ③ with the adjuster screw ④ and then re-tighten the counter nut (fig.24).

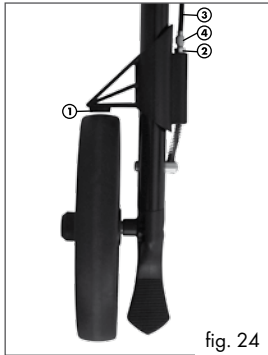Tyres

1. Heavily worn tyres, or tyres with visible damage (e.g., cracks), must be exchanged. Do not use the Rollator if this damage is apparent.
2. Check if all of the wheels are turning freely at least once a week.

Frame

Check the frame of the Rollator for cracks and damage on a regular basis.

Repairs

Never attempt to change any of the product characteristics and do not attempt to repair this article yourself! Obtain spare parts only from the manufacturer. Unforeseeable damage can occur from improper repair work!

Disposal

Please dispose of the product and packaging in an environmentally friendly manner, separating different materials for recycling as required! Dispose of this item through an authorised disposal company or through your local waste disposal facilities - ask your local authority for information if further guidance is required. Ensure that you comply with all regulations currently in force.

3 Year guarantee

This product was manufactured under constant quality controls. We guarantee that it is free from manufacturing and material defects. The guarantee claim is not valid for normal wear and tear and is invalid in case of misuse.

Should you however have cause for complaint, please contact the free phone service hotline, visit our website, or send us an e-mail. You will receive personal advice in all cases. During the guarantee period, you can arrange with our service staff that you send the faulty component(s) to the stated service address. We will immediately replace or repair the component, and send it back to you. Please retain the till receipt. Describe the nature of your complaint, and we will do our best to find a satisfactory solution as soon as possible.

IAN: 46438

Service address

Product hotline 00800-88080808 (Mon.-Fri. 10am – 4pm, free phone)

(GB) PSB LOGISTICS MANAGEMENT LTD • UNIT B3 SITE B
THE BOWMAN CENTRE ARNCOTT, OX25 1NZ OXON

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Ⓒ Atlantic Forwarding Ltd. • Mühlemattstr. 22, C.P. No. 131